
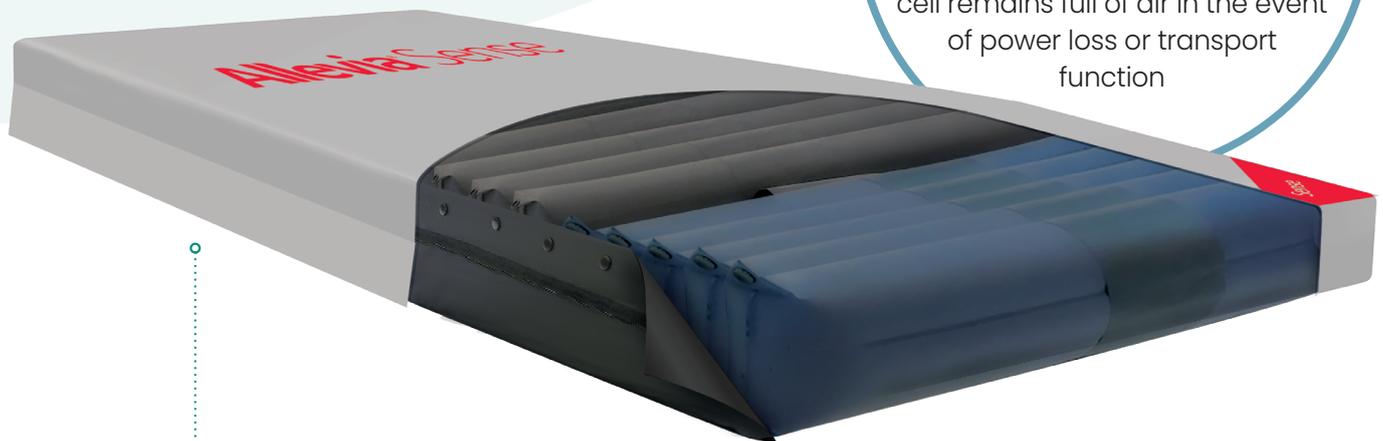


# Allevia<sup>®</sup>Sense

## Suitable for high to very high risk residents

This mattress uses a dynamic mattress system with alternating cells for pressure relief. The cells are 18 cm deep to provide deep immersion, and the heel slope is created by lower heel cells. The auto weight pump features an adjustable comfort setting with the option of constant low pressure therapy. This mattress also features an ultra rapid CPR valve.

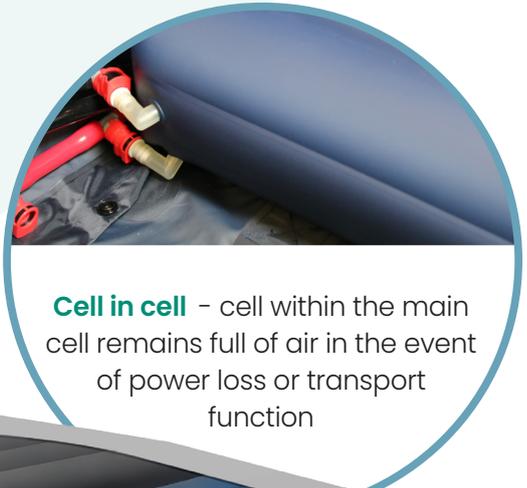

Size	2000 x 900 x 190 mm (79 x 35 x 7.5") 2000 x 1050 x 190 mm (79 x 41 x 7.5")
Weight capacity	249 kg (39 st)
Warranty	2 years
Cover	Four-way stretch PU
Vapour permeable	✓
Layered foams	n/a
Zip	✓
Water resistant	✓



**Auto weight /**  
pressure set pump  
unit