

# 600 Mattress

## Maximizing support as well as comfort

Pressure support surface with gel-infused open cell memory foam in head and torso sections to provide **ample support** while keeping the resident **comfortable**.

### Advanced pressure support, delivering comfort and relief

Superior pressure redistribution, delivered by enhanced support - means less cradling at higher levels of compression and eliminating static fatigue and shear.

Protects vulnerable heel area with 30-degree slope - redistributing pressure to the entire leg, away from the heel.

Firm side bolsters provide additional support for easy patient transfer.

### Measurements:

- Dimensions:  
36" x 80"  
36" x 84"  
42" x 80"  
42" x 84"
- Weight capacity: 525 lbs

### Materials & Certification:

- High-performing foams: Yes
- Cover: 4 way stretch nylon top cover reduces perspiration and provides a low shear friction surface, that is anti-microbial and fluid resistant
- Bottom cover: Triple laminate vinyl, fluid resistant, anti-microbial/anti-bacterial, fire retardant
- Cover color: Blue
- Vapor permeable cover: Yes
- Zipper: 360-degree zipper for easy removal, cleaning or replacement. Anti-Contamination Flap covers the zipper, reducing the possibility of fluid migration
- Waterproof: Yes
- Fire retardancy: Fire sleeve meets Federal Fire Code 16 CFR part 1633
- ACC600-0-FMI-100: Standard mattress
- ACC600-0-FM2-100: Wide mattress

### Warranty:

- Warranty: 5 years for foam, 2 years for cover

#### 1. Side Bolsters

Firm high density foam

#### 2. Base layer

Firm high density foam circulation cuts to prevent tearing

#### 3. Middle layer



Firm medium density foam

#### 4. Top layer

Gel infused memory foam

#### 5. Heel Slope

Soft foam

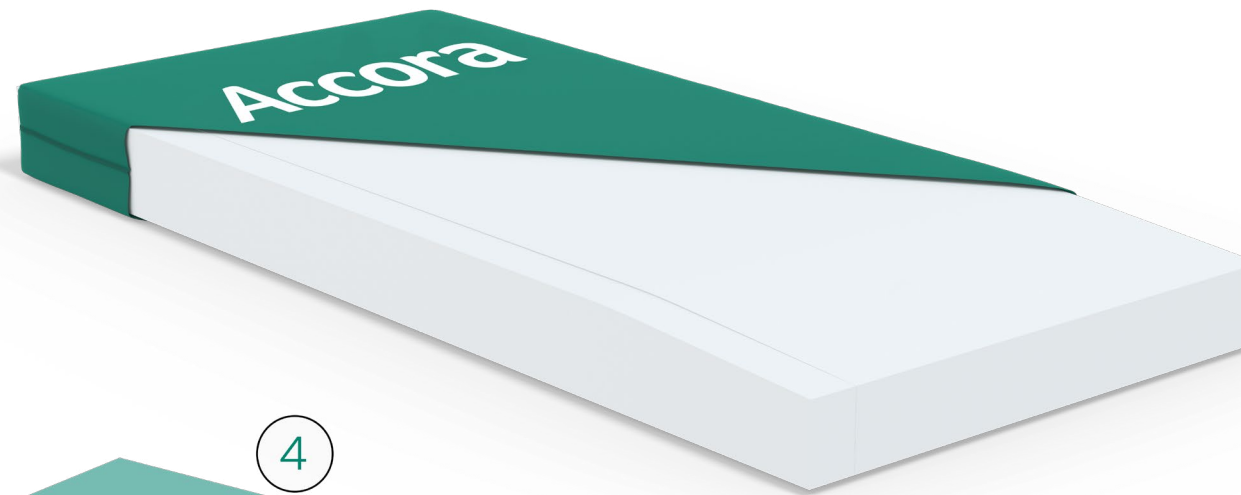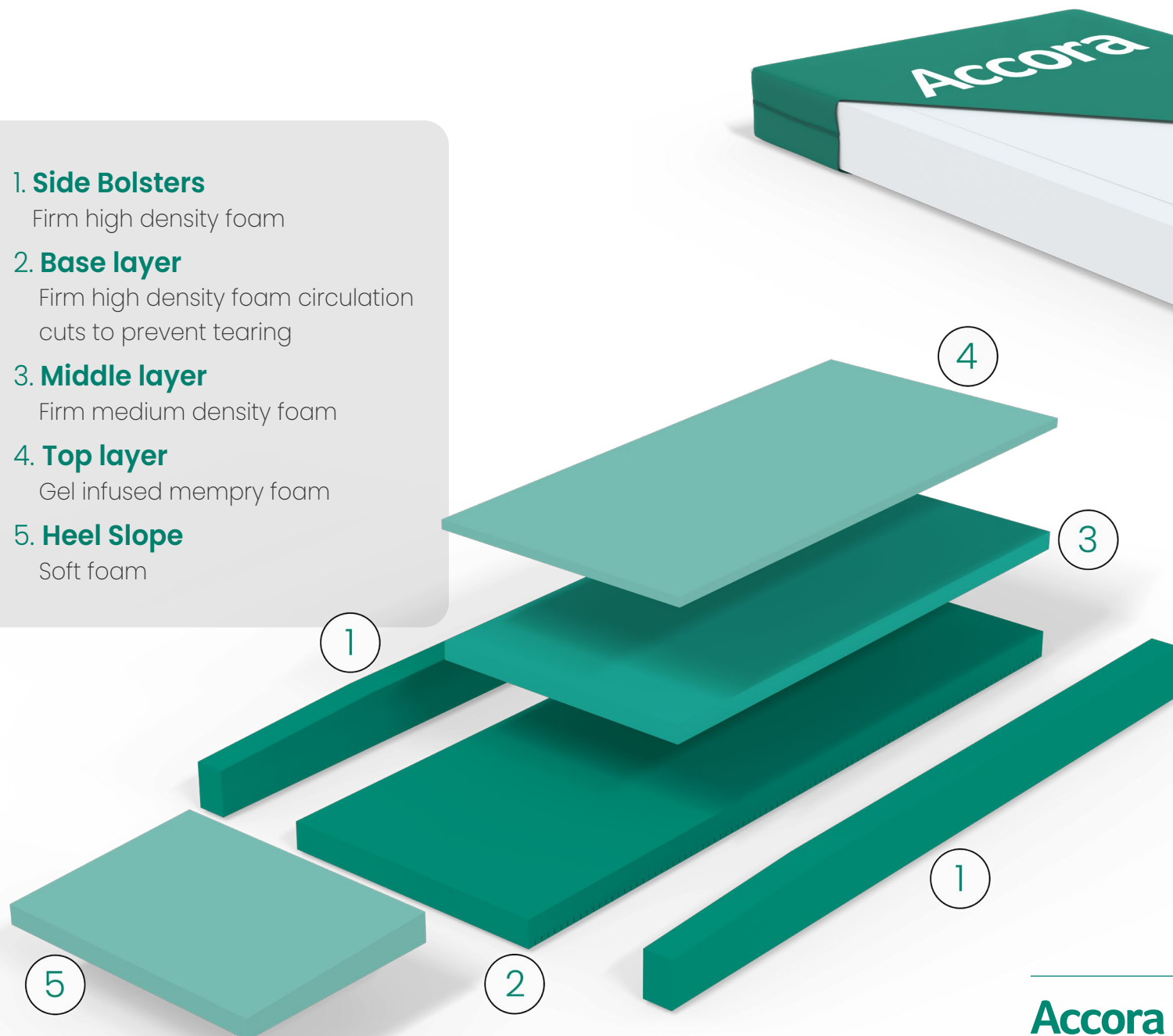

Provide **great clinical support** and **comfort** for your residents

Be compliant. Be **safe**. Be functional

Pioneers in FloorBed® Technology

Accora

6021 University Blvd, Suite 530, Ellicott City, MD 21043

T: 301 560 2400

E: [information@accora.care](mailto:information@accora.care)

W: [us.accora.care](http://us.accora.care)