


Allevia[®]


Pressure care solutions

Allevia® Classic

Reliable pressure redistribution
for daily comfort and protection.

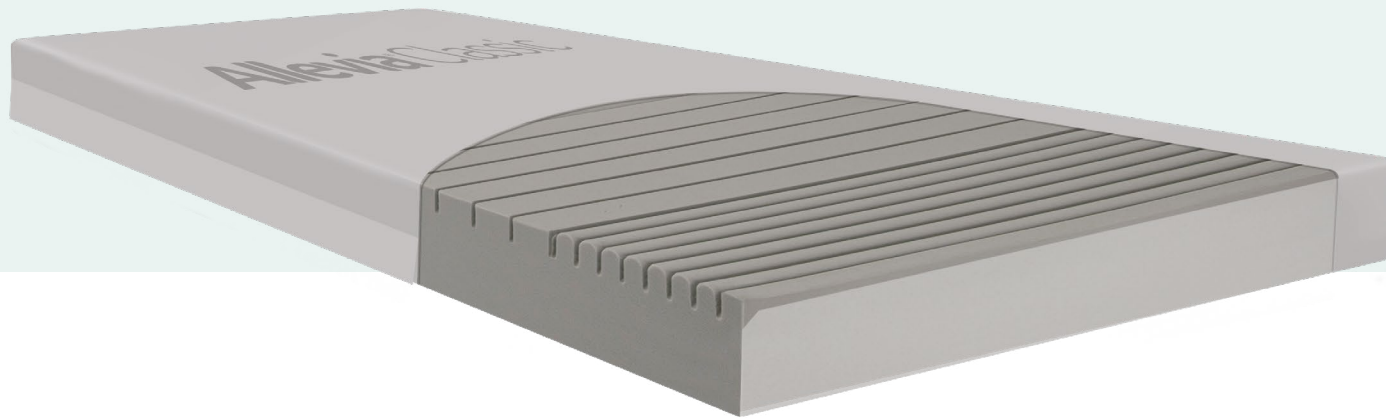

Size	36" x 80" / 42" x 80"
Mattress height	6 inches
Weight capacity	350 lbs
Cover	Top: Multi-directional stretch nylon, antimicrobial Bottom: Vinyl, anti-microbial/anti-bacterial
Zip	360-degree zipper for easy removal, cleaning or replacement
Fluid resistant	✓
Fire retardancy	Fire sleeve meets Federal Fire Code 16 CFR part 1633
Warranty	Foam: 4 years Cover: 2 years

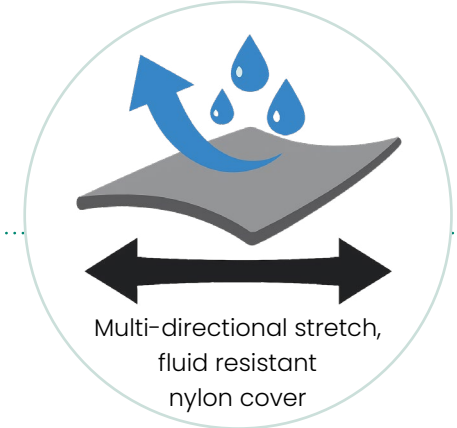


Allevia® Comfort


Provides consistent pressure distribution to help reduce the risk of injury.

The die-cut castellations in the high-density foam allow the mattress to contour closely to the body, promoting consistent pressure distribution helping to relieve high pressure points.

Firm perimeter promotes safer ingress and egress, aiding mobility and helping mitigate fall risk.

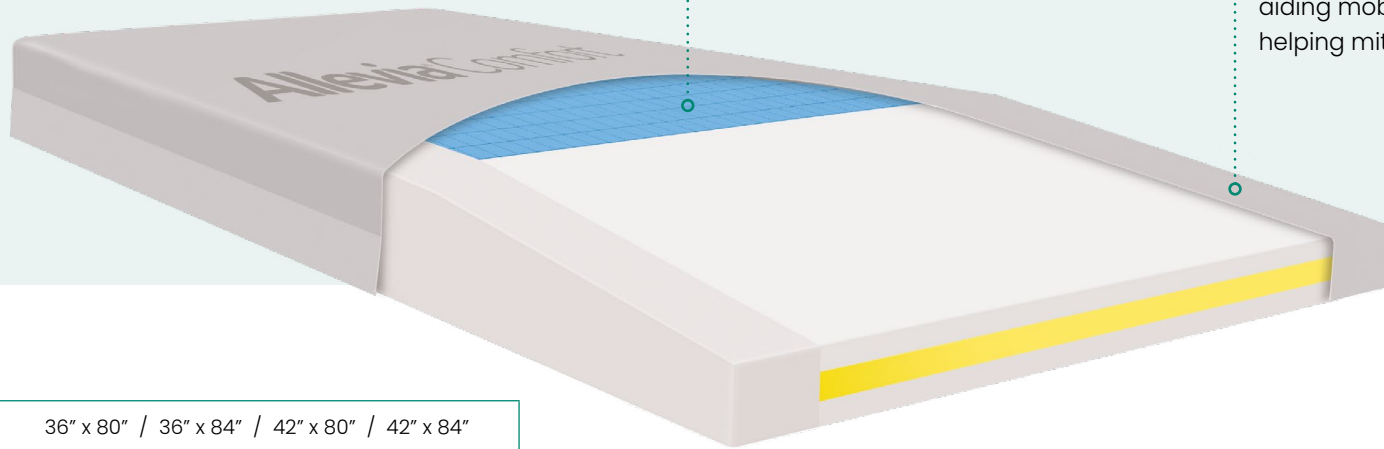


Size	36" x 80" / 36" x 84" / 42" x 80" / 42" x 84"
Mattress height	6 inches
Weight capacity	400 lbs
Cover	Top: Multi-directional stretch nylon, antimicrobial Bottom: Vinyl, anti-microbial/anti-bacterial
Zip	360-degree zipper for easy removal, cleaning or replacement
Heel slope	✓
Firm edges	✓
Fluid resistant	✓
Fire retardancy	Fire sleeve meets Federal Fire Code 16 CFR part 1633
Warranty	Foam: 5 years Cover: 2 years

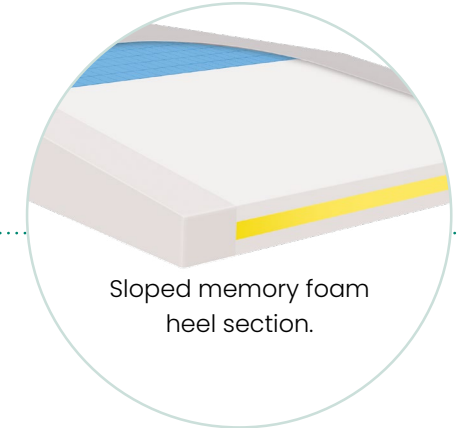


Allevia® Prima


Rapid-cooling, gel-infused foam that helps preserve skin integrity.

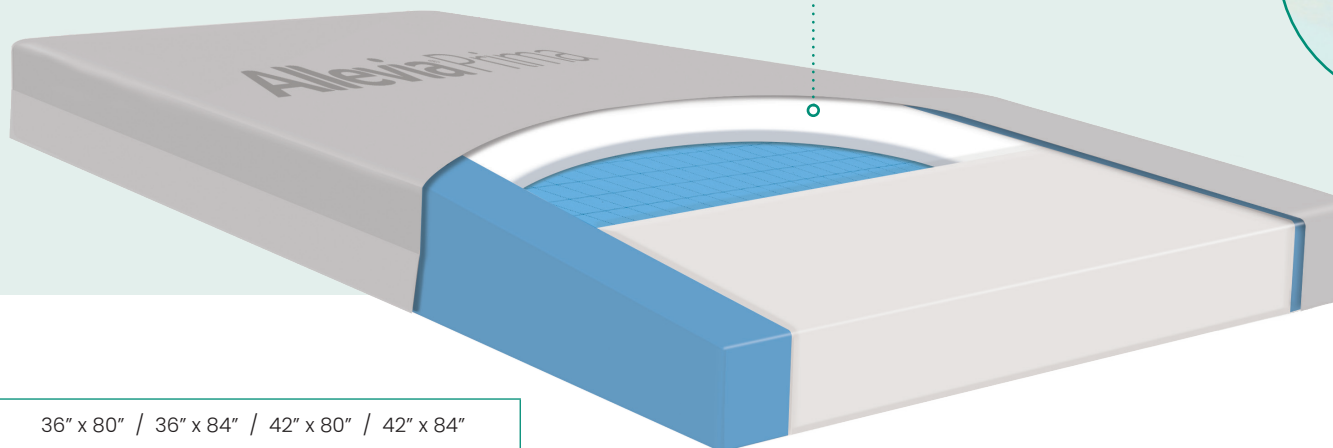


Gel-infused memory foam top layer draws heat away from the body, reduces moisture buildup, and helps maintain skin integrity.

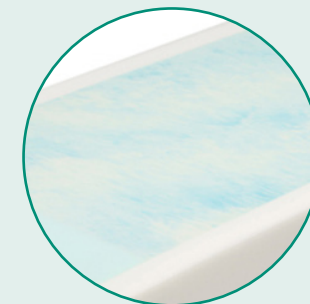

Size	36" x 80" / 36" x 84" / 42" x 80" / 42" x 84"
Mattress height	6 inches
Weight capacity	500 lbs
Heel slope	✓
Firm edges	✓
Cover	Top: Multi-directional stretch nylon, antimicrobial Bottom: Vinyl, anti-microbial/anti-bacterial
Zip	360-degree zipper for easy removal, cleaning or replacement
Fluid resistant	✓
Fire retardancy	Fire sleeve meets Federal Fire Code 16 CFR part 1633
Warranty	Foam: 10 years Cover: 2 years

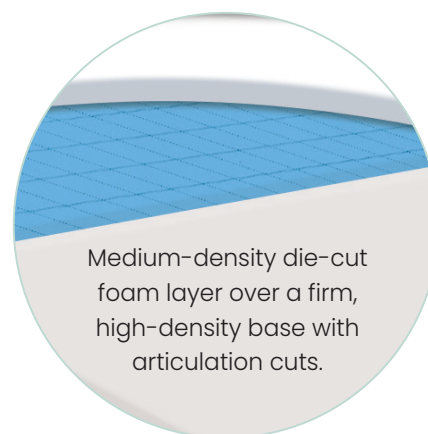

Medium-density die-cut foam layer over a firm, high-density base with articulation cuts.

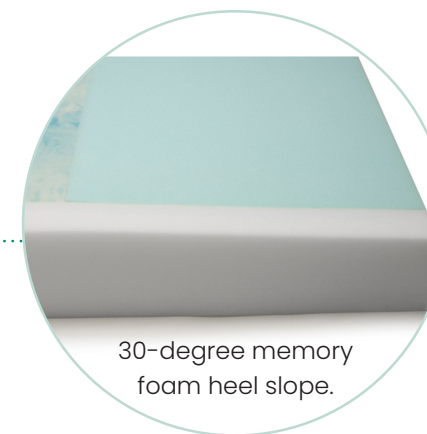


30-degree memory foam heel slope.

Allevia® Duo


Self-adjusting technology that responds to movement.

Delivers non-powered, self-adjusting zoned pressure redistribution, and converts to a dynamic powered air mattress with the use of the optional pump.

Indicated for the prevention and treatment of any and all stage tissue injury

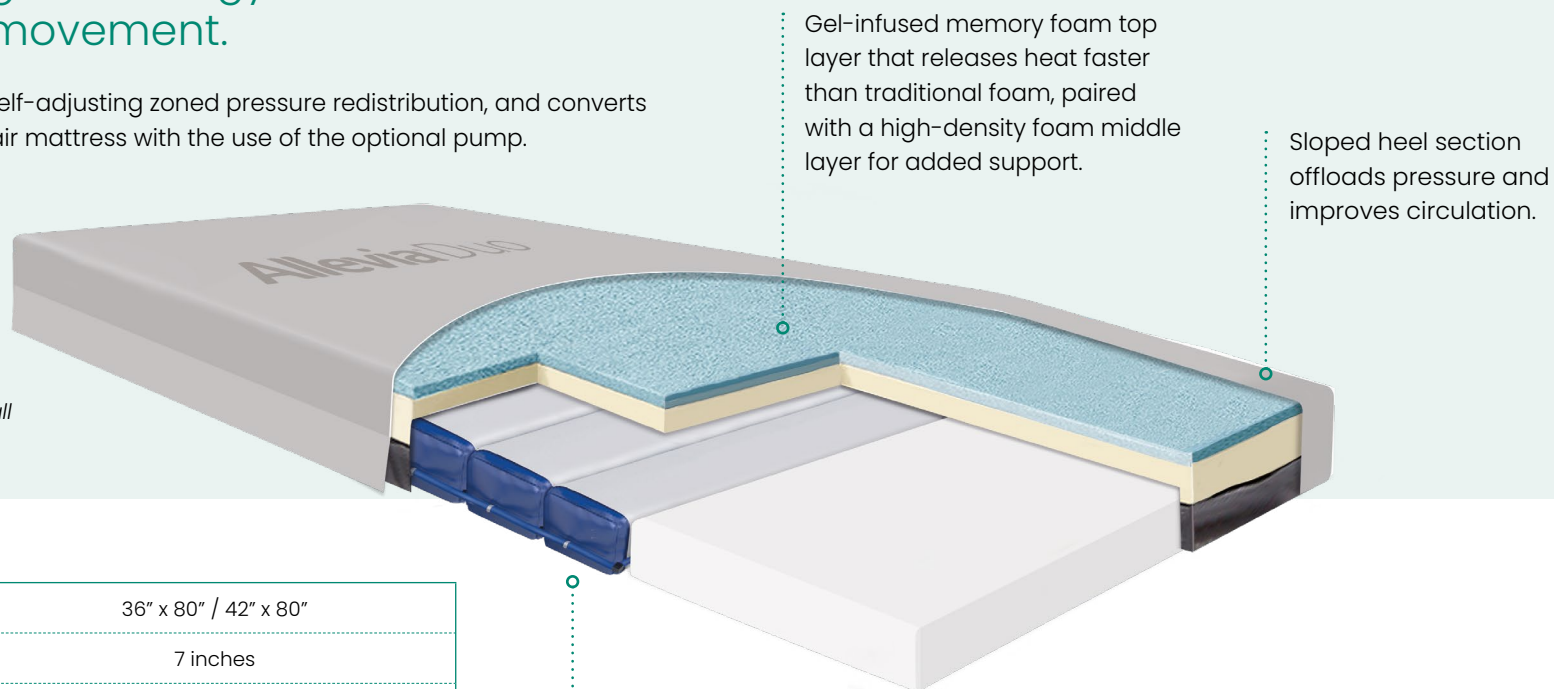

Size	36" x 80" / 42" x 80"
Mattress height	7 inches
Weight capacity	36" - 500 lbs 42" - 600 lbs
Heel slope	✓
Cycle time	10 minutes
Cover	4-way stretch nylon, antimicrobial Bottom: Vinyl, anti-microbial/anti-bacterial
Zip	360-degree zipper for easy removal, cleaning or replacement
Fluid resistant	✓
Fire retardancy	Fire sleeve meets Federal Fire Code 16 CFR part 1633
Warranty	5 years for mattress 18 months for pump


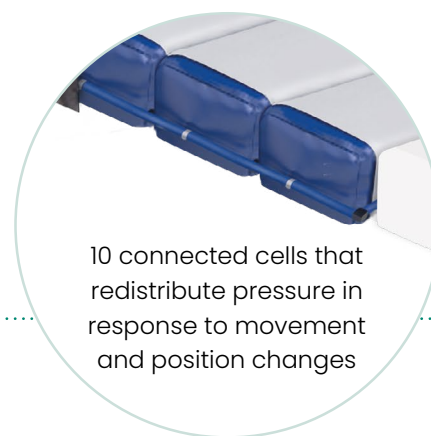


Allevia[®]Sense


Individualized pressure care with adjustable air cells and Auto-Sense™ weight detection.

7-inch cells offer deep immersion and can be repositioned or deflated to relieve pressure on areas such as the kyphotic spine.

Auto-sense and pressure set pump unit.

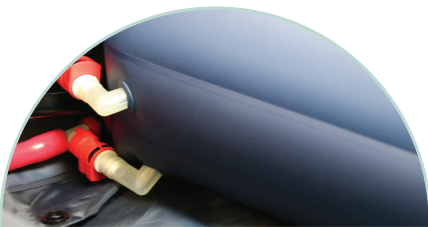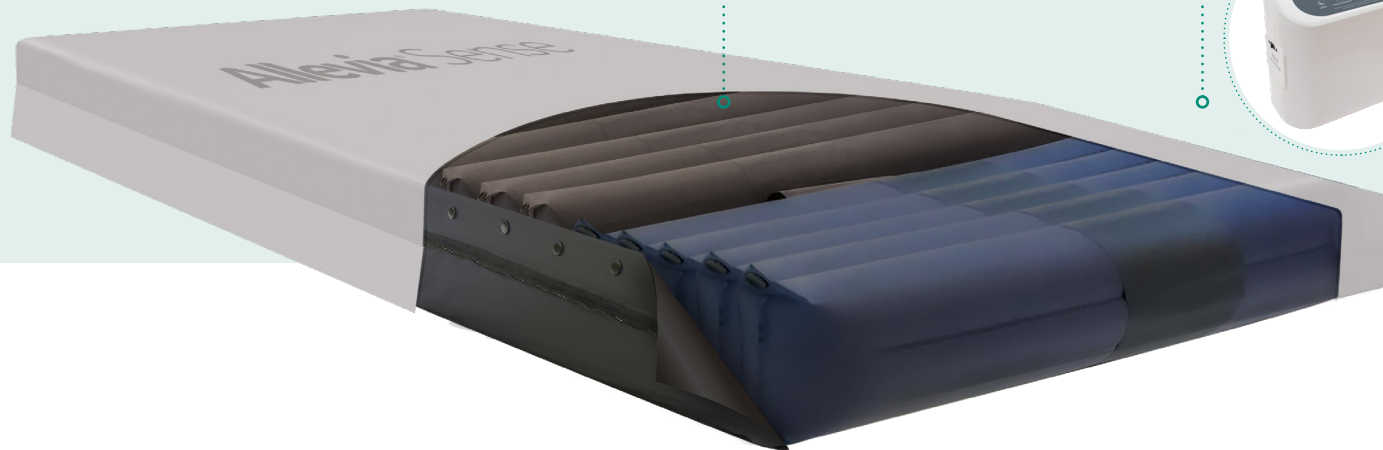

Cell in cell – cell within the main cell remains full of air in the event of power loss or transport function.

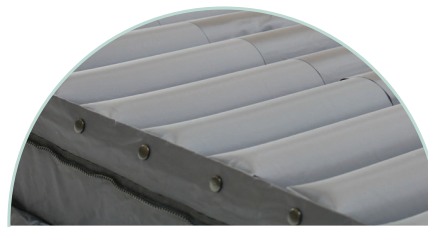

Oval top cells for even and complete pressure redistribution.




Ultra-rapid CPR valve.



Heel slope – 4.5-inch cells create a heel slope, with the option to adjust for further offloading on the heels.

Auto-sense pump features automatic weight monitors – no need for manual weight input. Option for manual override to adjust for patient comfort.

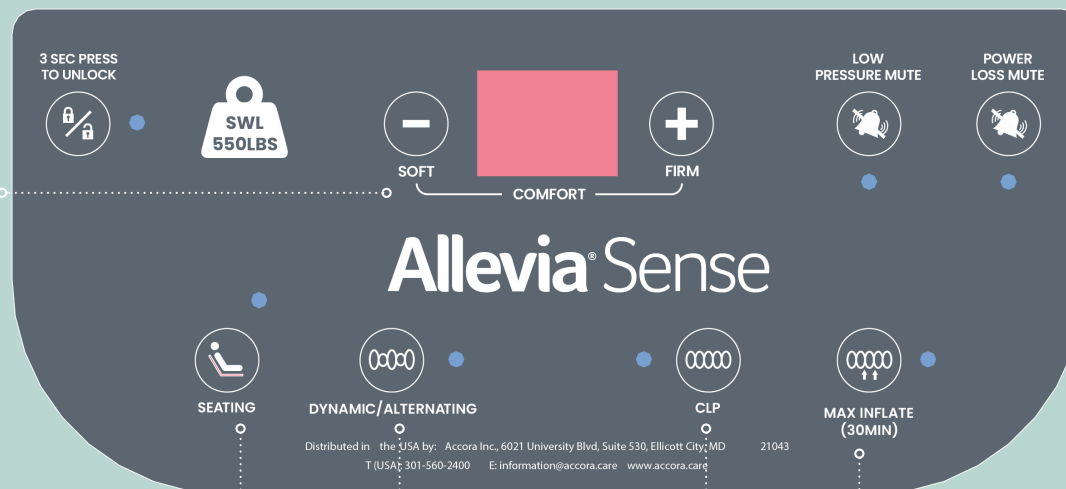

Seat inflate option
– adjust pressure for when resident is in upright position

Dynamic
1 in 2 cycle alternating therapy

Constant Low Pressure
therapy action

Autofirm mode gives static maximum inflation to allow nursing and care procedures and auto turns off after 30 minutes

Sizes	36" x 80" / 42" x 80"
Weight capacity	550 lbs
Cycle time	10 minutes
Cell type	2 in 1
Cover	4-way stretch fluid resistant PU cover
Zip	360-degree zipper for easy removal, cleaning or replacement
Fluid resistant	✓
Fire retardancy	Mattress cover is fire-retardant
Warranty	2 years for mattress and pump

Accora

Accora Inc,
6021 University Blvd,
Suite 530,
Ellicott City,
MD 21043

T: 301-560-2400
E: information@accora.care
W: us.accora.care

Accora

Accora Inc,
6021 University Blvd,
Suite 530,
Ellicott City,
MD 21043

T: 301-560-2400

E: information@accora.care

W: us.accora.care