

# Accora BodyMove™ Technology

Designed to reduce pressure and shear forces  
by moving with the body — not against it.

Available on  
the Empresa  
FloorBed




## What is BodyMove technology?


During profiling, **BodyMove** shifts the backrest upward and backward, and moves the thigh platform toward the foot of the bed.

This multi-directional movement:

- Reduces shear forces
- Supports patient comfort
- Supports patient handling practices and may reduce staff strain

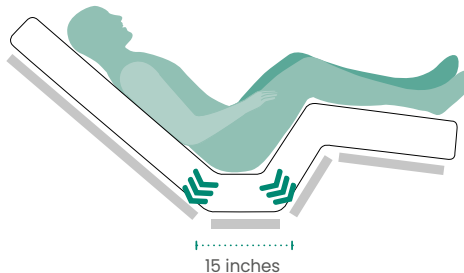


**Empresa** with **BodyMove** regression reduces friction and shear, and aids pressure distribution

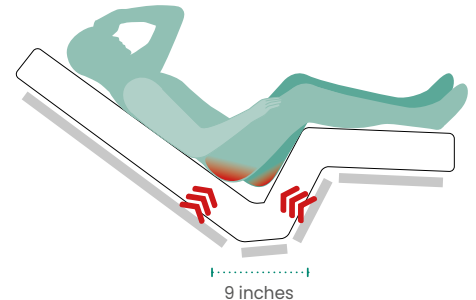


Standard 4 section sleep decks with no regression increases friction and shear

## Study: Empresa with BodyMove vs standard bed

A study assessed the impact of BodyMove technology on sacral pressure by comparing the Empresa bed with BodyMove to a standard nursing bed without it.

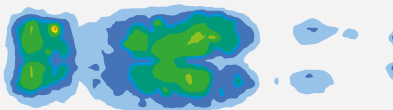

Pressure map for Empresa bed  
with BodyMove regression

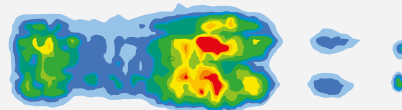

Pressure map for standard  
nursing bed with no regression

### The result:

Average pressure over  
the sacral area was

reduced by  
**23%**

when using the Empresa  
bed compared with the  
standard nursing bed.

**Accora. Pressure mapping white paper:** Assessing the effect of Empresa® bed BodyMove technology on sacral interface pressure compared with a standard nursing bed. Available at [us.accora.care](http://us.accora.care)